

Siemens Healthcare Diagnostics GmbH, SHS EMEA CEET QT, Siemensstrasse 90,
1210 Vienna

Name M.A. Roland Ertl
 Department SHS EMEA CEET QT
 Telephone +43 51707-38274
 Mobile +43 (664) 8011738274
 E-mail roland.re.ertl@siemens-healthineers.com
 Date January 14, 2020
 Document Ref# ASW20-02.A.OUS

Urgent Field Safety Notice:

Atellica® IM 1300 Analyzer / Atellica® IM 1600 Analyzer

Humidity Packs Incorrectly Identified as Expired (Lot 0010)

Dear Sirs,

Our records indicate that your facility may have received one or more or a combination of the following products:

Table 1. Atellica® Solution Affected Product(s):

Product	Siemens Material Number (SMN)
Atellica IM Humidity Pack (Qty 5)	11313505
Atellica IM Humidity Pack (Qty 1)	11313496

Reason for Urgent Field Safety Notice

The purpose of this communication is to inform you of an ongoing investigation involving the products listed in Table1 above and provide instructions on actions that your laboratory must take.

Siemens Healthcare Diagnostics has received customer complaints stating that Atellica IM Humidity Packs Lot 0010 have expired, causing the analyzer to eject the Humidity Packs and stop processing samples.

Atellica IM Humidity Packs do not have a shelf life (lot) expiration, therefore, the system is identifying these packs as expired in error. Preliminary investigation indicates this error is attributed to the software.

It is important to note that Humidity Packs do have an onboard stability (OBS) of 180 days that is monitored by the system software. Once loaded onto an analyzer, the OBS of each pack is monitored.

Siemens is actively working to investigate the root cause and customers will be notified when additional information is available. This behavior will be corrected in a new version of software that will be available shortly.

Risk to Health

The potential exists for an apparent delay in testing when this issue occurs. Siemens is not recommending a review of previously generated results as the accuracy of results is not affected by this issue.

Actions to be Taken by the Customer

The error can be resolved by manually entering the Humidity Pack barcode information with the following steps.

1. Cover the existing Humidity Pack 2D barcode with permanent marker so it cannot be read by the system. Do not remove the original barcode. (top view)



2. On the Command bar, select **Inventory**.

Note: Steps 2 through 8 correspond to the screen shot on page 3 below.

3. Select the **Reagent Loader** tab.

- a. Select the **IM module** where the packs are to be loaded.

4. Select **Manual Entry**.

5. Select pack type as **Ancillary**.

6. Enter the following barcode number.

- a. 78710600PPPPPP

- b. PPPPPP is a unique six-digit pack number. Use the unique identifier on the Humidity Pack being manually entered.


- c. The unique identifier on the Humidity Pack label is depicted below. In this example the unique identifier is 034831. (side view)



7. Select Verify.

8. Confirm Assay equals Humidity and the date next to the hourglass symbol is 2020-02-29.

9. Select OK.
10. Load the Humidity Pack in position 5 of the reagent tray.
11. Repeat this procedure for all required Humidity Packs. Each analyzer will require, 7 to 10 packs based on analyzer environmental requirements. When the full set of required Humidity Packs are loaded, the analyzer will resume normal operation.

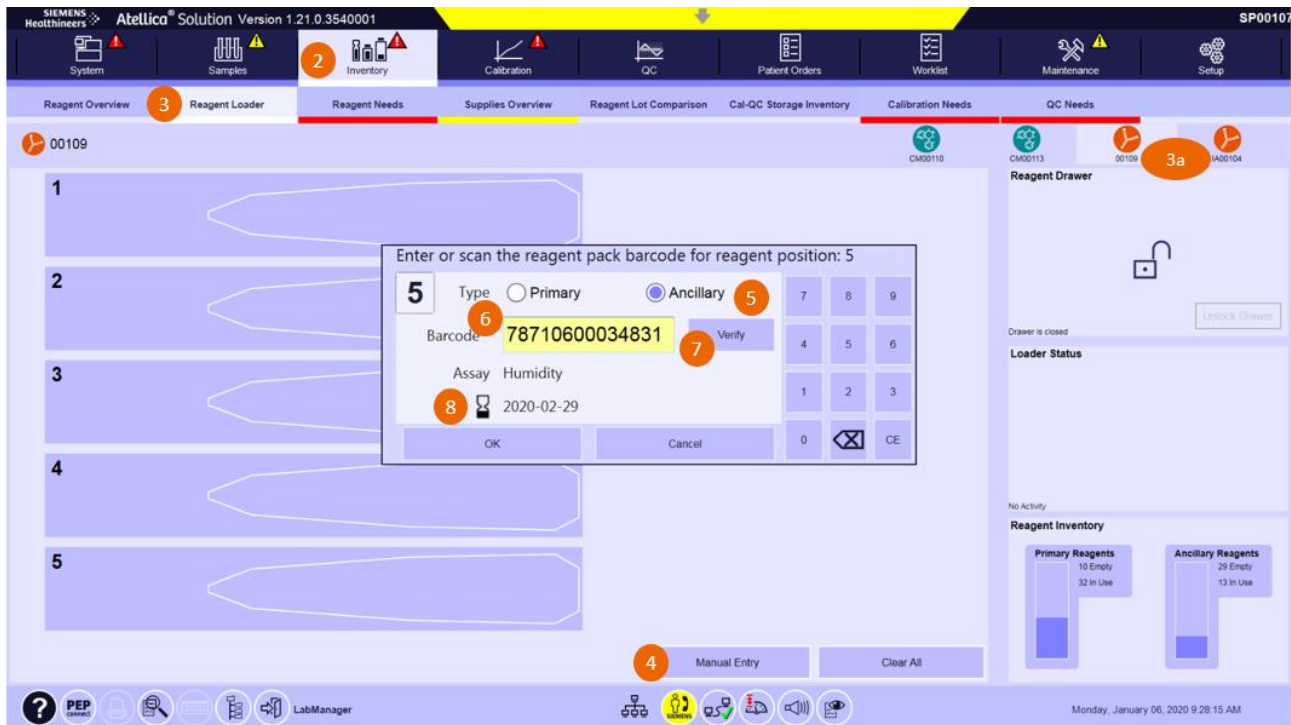

Note: Anytime the analyzer must repeat the “Startup” routine, the humidity packs with covered barcodes will be identified as “unknown.” This will require the above actions to be completed again. The analyzer will perform the “Startup” routine after a reboot, mechanics off/on and/or opening of the front or rear top covers.

- Please review this letter with your Medical Director.
- Complete and return the Field Correction Effectiveness Check Form attached to this letter within 30 days.
- If you have received any complaints of illness or adverse events associated with the products listed in Table 1, immediately contact your local Siemens Customer Care Center or your local Technical Support representative.

Letter of January 14, 2020
to



Please retain this letter with your laboratory records and forward this letter to those who may have received this product.

We apologize for the inconvenience this situation may cause. If you have any questions, please contact your Siemens Customer Care Center or your local technical support provider.

Sincerely yours,

Siemens Healthcare Diagnostics GmbH

A handwritten signature in black ink, appearing to read "Franz Schwarz".

i.V. Dipl. Ing. Franz Schwarz
Quality Management CEE

A handwritten signature in blue ink, appearing to read "Gernot Osterer".

i.A. Ing. Gernot Osterer
Product Manager Austria & SEE